



App Quick Start Guide


Discovery RA App

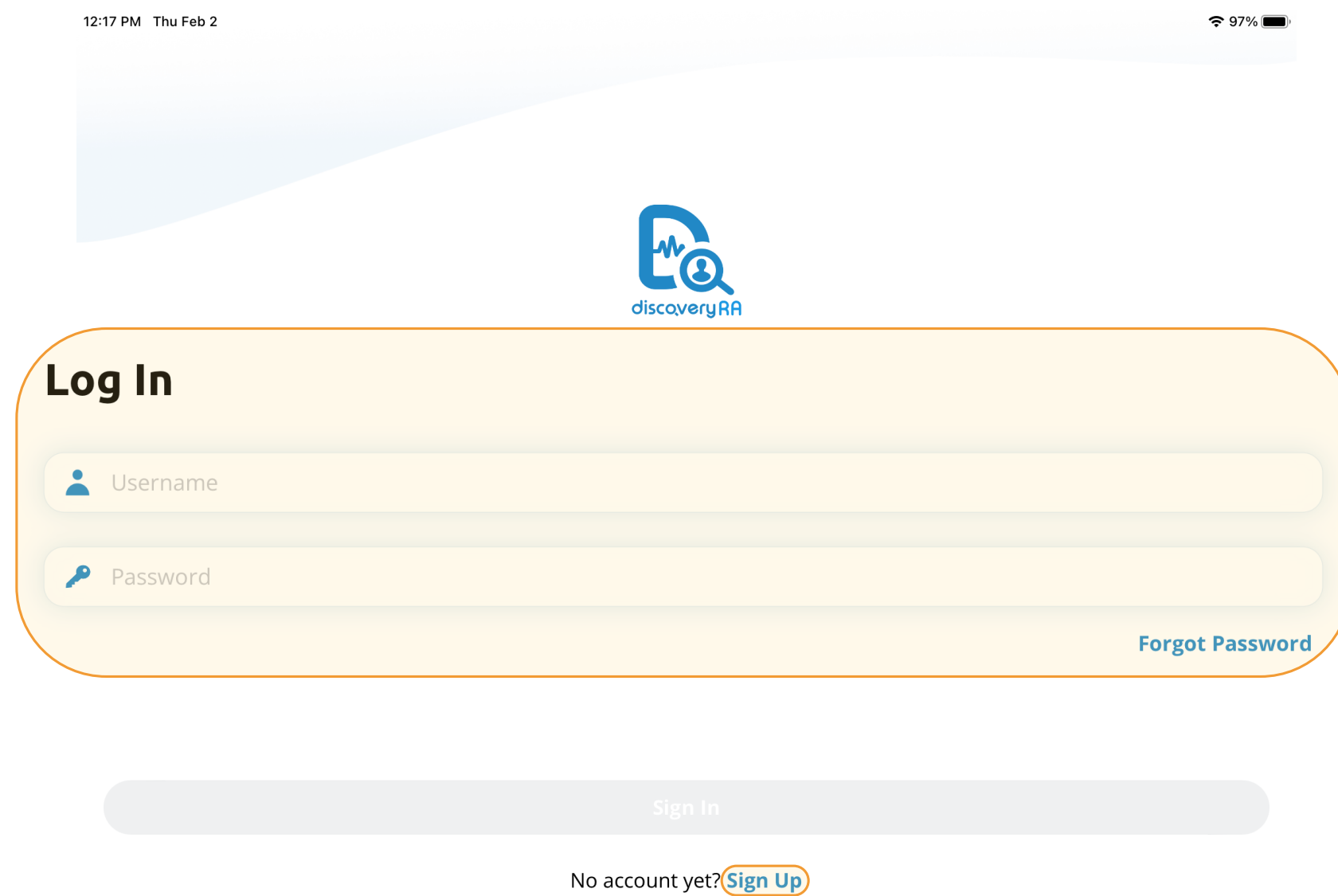


Log in to the Discovery RA app by entering a username and password. You can create an account if you do not have one, only an email is required.

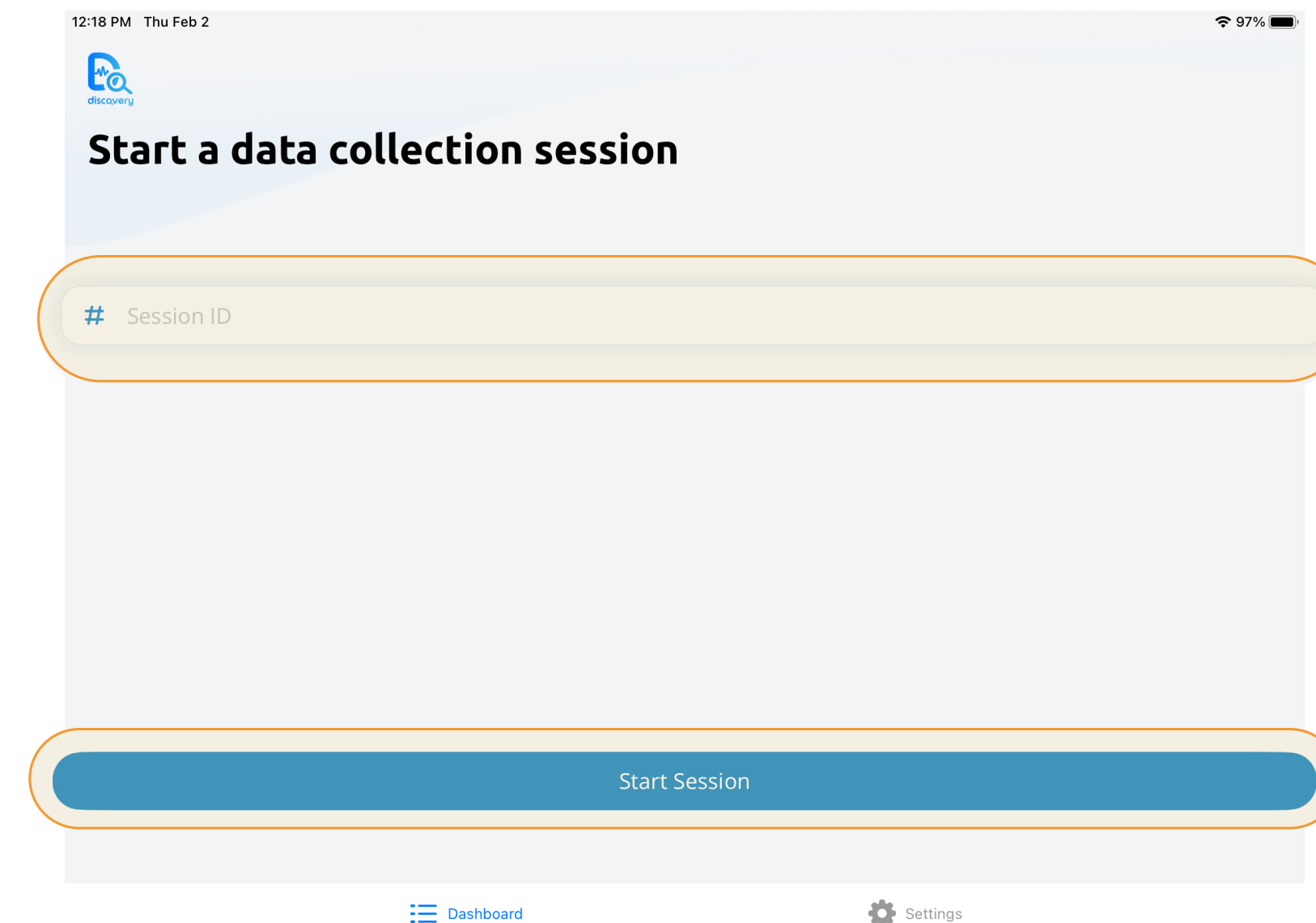

Type in a number to represent the data collection Session ID to be used. Then click “Start Session” to begin.

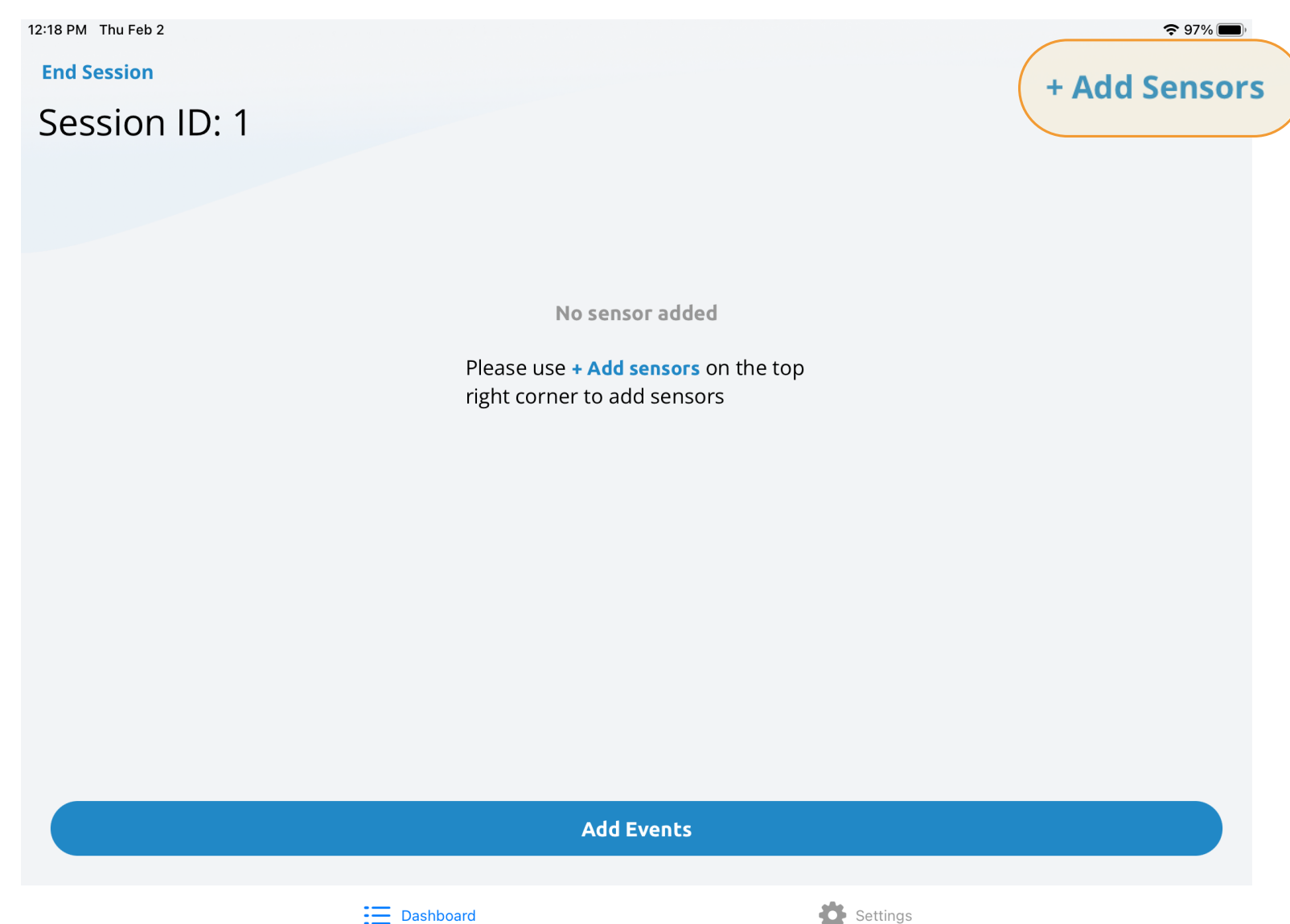


Click “+ Add Sensors” on the top right corner to display the list of advertising devices.

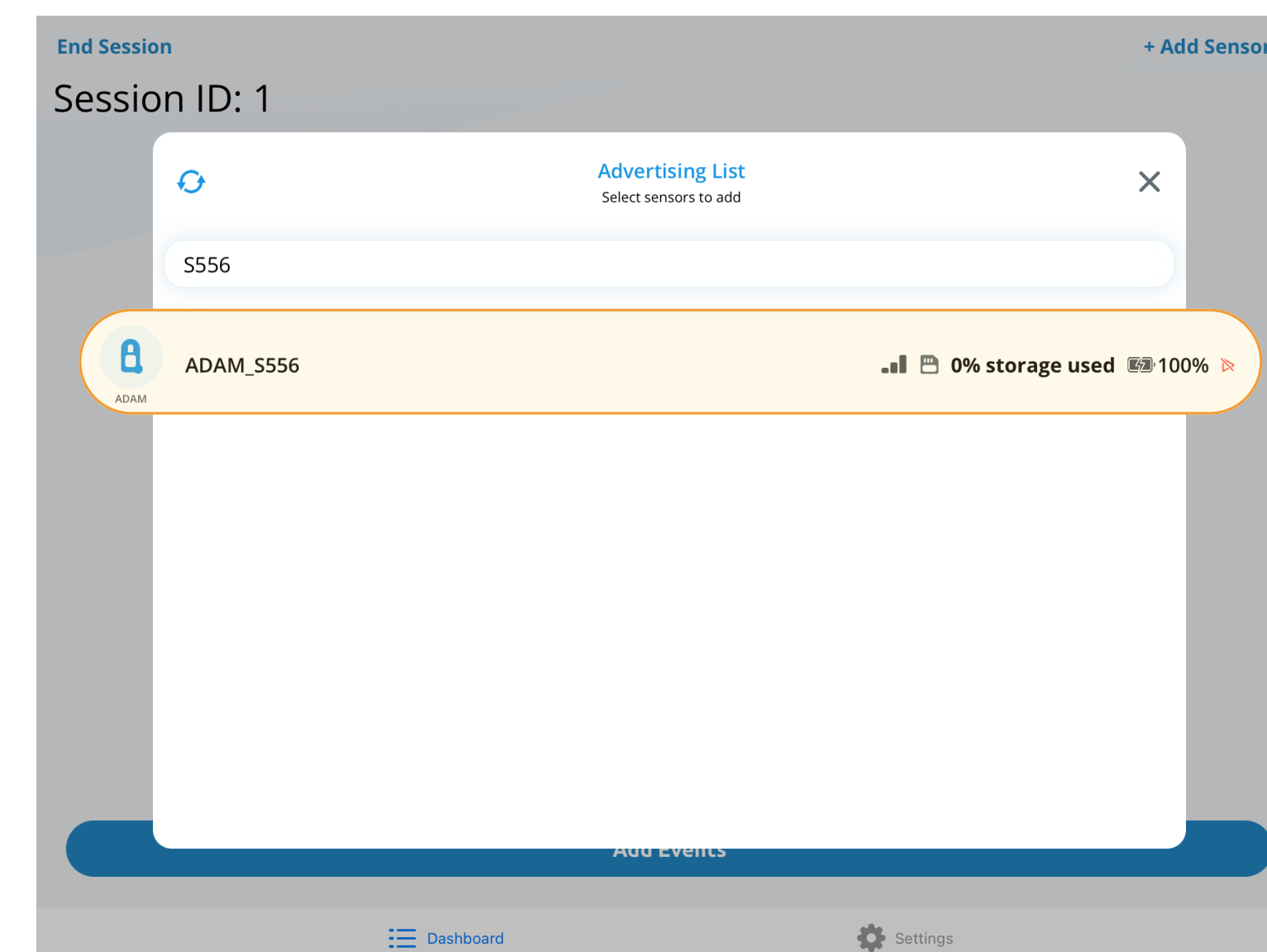

Click on the sensor and a window will appear, prompting confirmation to connect to the sensor. If prompted to pair the device via Bluetooth, select “Pair”.

Discovery RA App

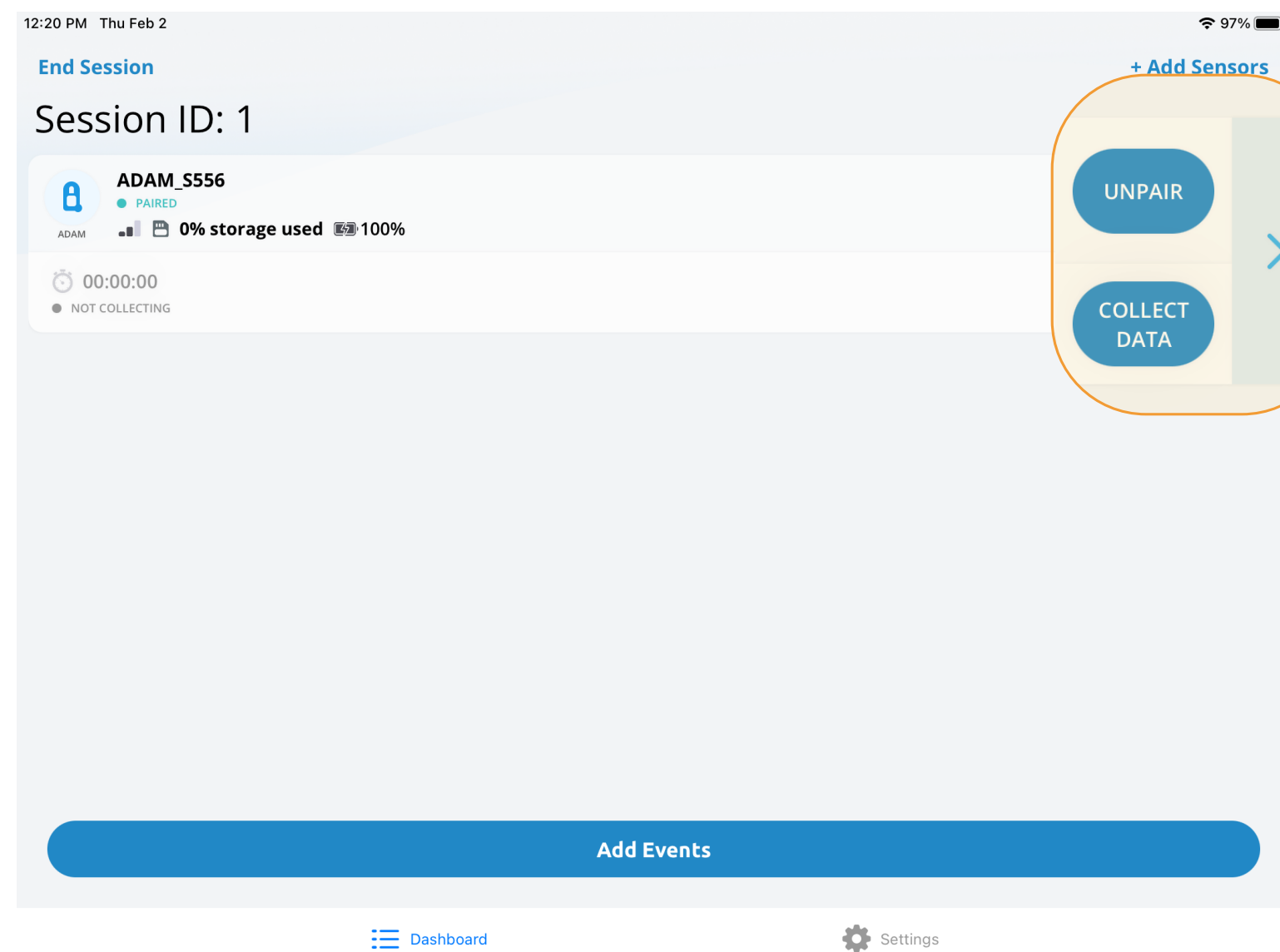


From the dashboard, click the right arrow on the sensor row to display a window with the sensor's information. Unpair the sensor by clicking "Unpair". Click "Collect Data" to begin recording data. Click "Stop" to end data collection.



Click on "Stream" to view live streamed data. The type of data displayed will depend on the capabilities of the sensor. Data collection can also be initiated and stopped here by clicking "Collect Data" and then "Stop".



Click "Download" to save all recorded data off of the sensor and onto the mobile device and then to the cloud.

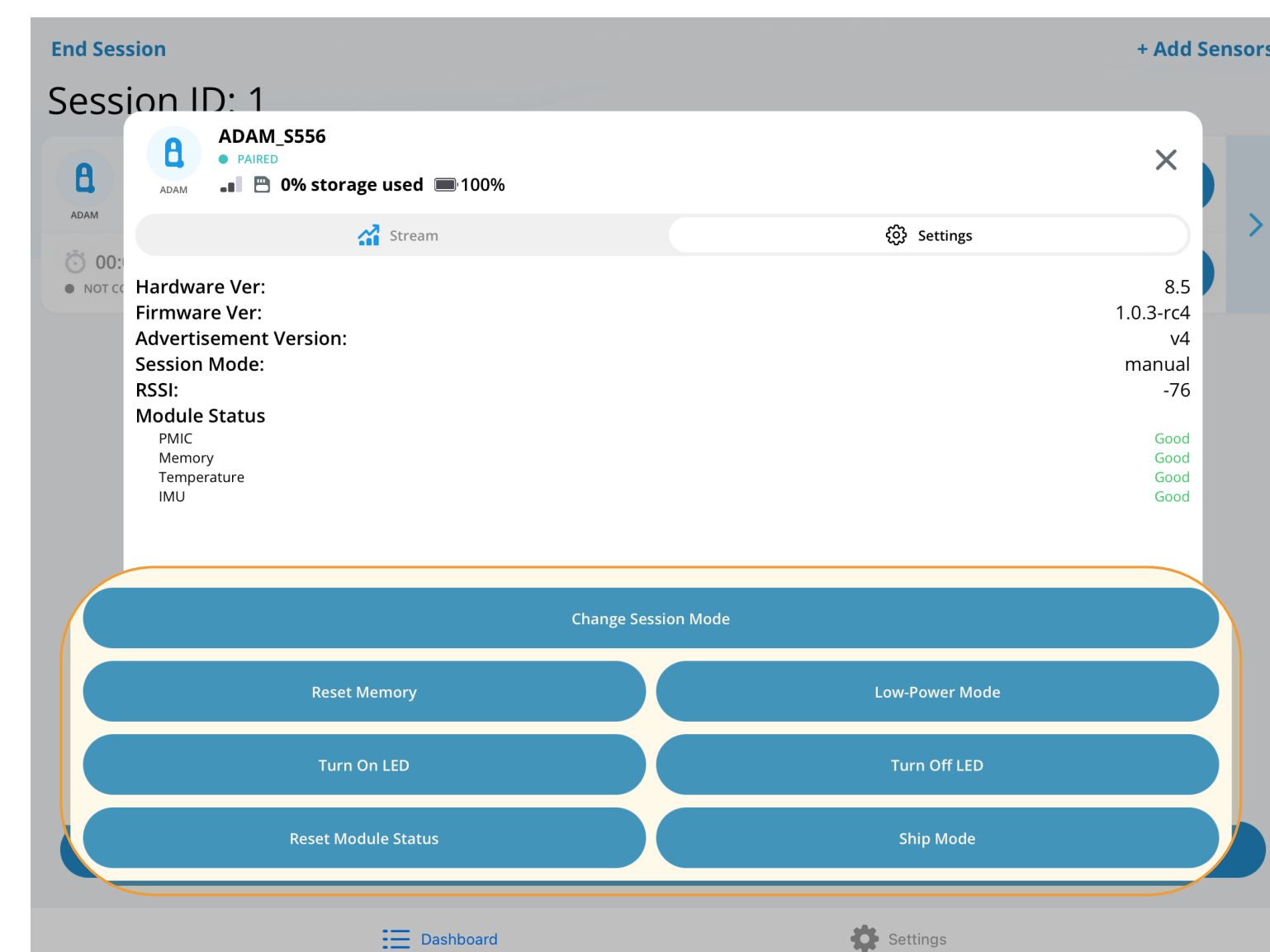

Click on the "Options" section for sensor options. Click "Reset Memory" to delete any data from the sensor. Click "Ship Mode" to turn off the sensor.