


Quick Start Guide

Daily routine



Complete the daily questionnaire if required

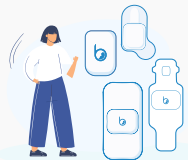

Wear the sensor



Collect the data




Sync and charge the sensor


How to setup the sensor and iPad

SENSOR SET UP

1 Sanitize the sensors and wireless charger.

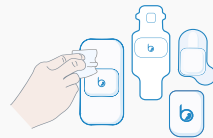


Sensors

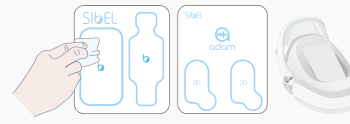

Wireless Charger

2 Place the sensors on the charger to wake up the sensors.

Anne, Anne Limb, Adam

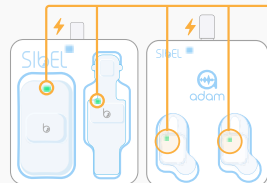

⌚ ~5 hours

It takes 5 hours to fully charge the sensor from 0%.

- Solid Green LED
Sensor is fully charged.
- Solid Red LED
Sensor failure.
- Flashing Blue LED
Sensor not connected to Discovery app after power-on.
Data available for download and on charger.
- Flashing Green LED
Sensor is charging, but not fully charged.
- Flashing Red LED
Sensor lead-off, low battery, or low memory.

Aria




- Solid Green LED
Sensor is fully charged.
- Solid Orange LED
Sensor is charging, but not fully charged.
- Solid Red LED
Charger failure.
- Solid White LED
Charger is functional and ready to use, but the sensor is not on the charger.
- LED Off
Charger is not plugged in.

IPAD SET UP

1 Make sure iPad is charged and connected to the internet.




2 Move the sensors close to the iPad to pair the sensors. (They have to be within 2 meters of each other.)

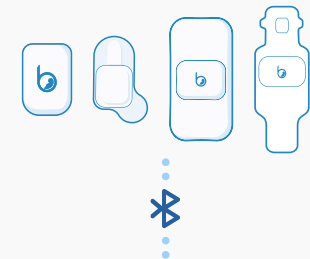

App Quick Start Guide

App manual

! Troubleshooting

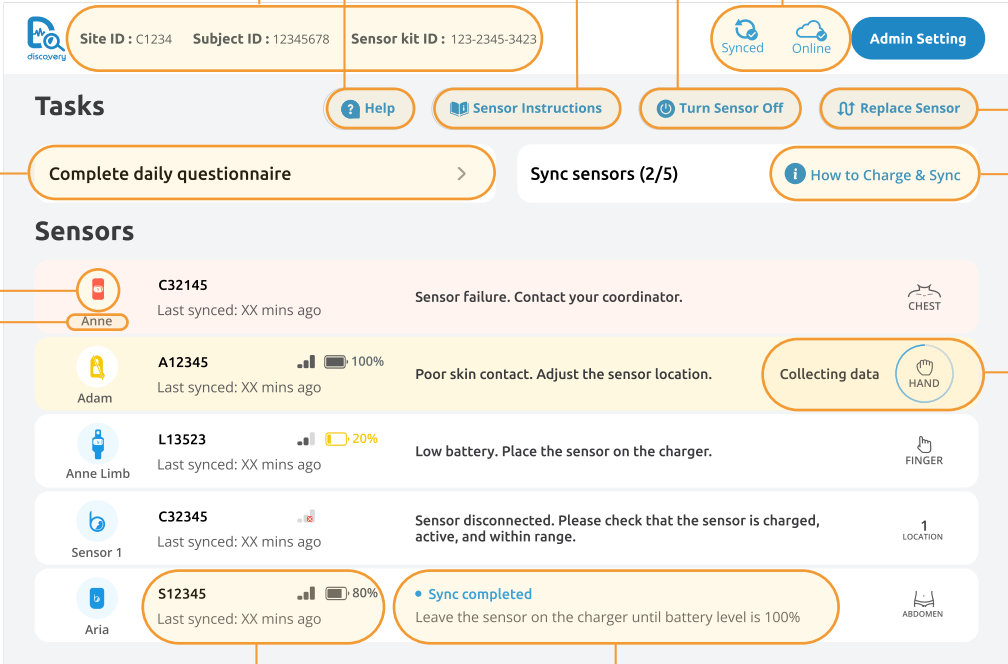

Help Center Button

Turn Sensor Off Button

Questionnaire Button

Kit Info

Site ID : C1234 Subject ID : 12345678 Sensor kit ID : 123-2345-3423

Sync/Online Status

Synced Online Admin Setting

Tasks

Help Sensor Instructions Turn Sensor Off Replace Sensor

Complete daily questionnaire > Sync sensors (2/5) How to Charge & Sync

Sensors

Sensor Status

Sensor Type

Replace Sensor Button

Charging & Syncing Instruction Button

Collecting Status / Body Location

Sensor Info

Syncing Status

Sensor ID	Subject	Status	Message	Body Location
C32145	Anne	Sensor failure	Sensor failure. Contact your coordinator.	CHEST
A12345	Adam	Poor skin contact	Poor skin contact. Adjust the sensor location.	HAND
L13523	Anne Limb	Low battery	Low battery. Place the sensor on the charger.	FINGER
C32345	Sensor 1	Sensor disconnected	Sensor disconnected. Please check that the sensor is charged, active, and within range.	1 LOCATION
S12345	Aria	Sync completed	Leave the sensor on the charger until battery level is 100%	ABDOMEN



Offline

Check the internet connection.



Unsynced

Check the internet connection.



Poor Skin Contact

Adjust the location or replace the adhesive.



Sensor Failure

Contact study coordinator.



Sensor Out of Range

Move tablet closer to the sensor. The sensor will automatically reconnect after the signal is detected.



Low Battery

Charge the sensor as soon as possible. The sensor will automatically disconnect when it is out of battery.